KOOGEEK


Instruction Manual

T1

———— Infrared Ear Thermometer ————

PACKAGE INCLUDES

1. Len
2. Probe
3. Battery cover
4. Temperature display mode: ear or forehead
5. Screen
6. Measurement and Memory Function button
7. On/off and Ear or forehead temperature Mode button
8. Protection/vertical storage cap.

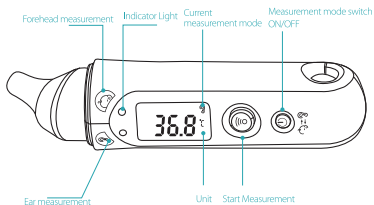UNIT (C / F) SWITCH

The product has two temperature units: Celsius $^{\circ}\text{C}$ and Fahrenheit $^{\circ}\text{F}$

Switch Methods: When the product is off, please press and hold "power" button until the display shows " $^{\circ}\text{C}$ " and " $^{\circ}\text{F}$ " alternately. Then release the button when the display is showing the unit you want.

MEASUREMENT MODES SWITCH

When the product is powered on, press "power" to switch between ear temperature mode and forehead temperature mode.

EAR MEASUREMENT

1. Remove the protection cap.
2. Press "power" button for 2 seconds to power on. The default mode is ear temperature mode. The indicator light of ear temperature mode will flash. When Smart Thermometer is ready for measurement, the indicator will stop flashing.
3. Ear measurement: Gently pull the ear upwards and backwards to straighten the ear canal, fix your head, and insert the probe at eardrum level. Press measurement button and finish the measurement when you hear notification sound.
4. Measurement results will be displayed on the screen and recorded automatically.

Ways for straightening ear canal:

1. Babies younger than 1 year old should pull the ear at rear direction.
2. Children and adults elder than 1 year old should pull the ear at outward and upper direction.

FOREHEAD TEMPERATURE MEASUREMENT

1. Remove the protection cap.
2. Press "power" button to power on, and continue to hold the power button to select the forehead mode. The indicator light of forehead temperature mode will flash. When Smart Thermometer is ready for measurement, the indicator will stop flashing.
3. Forehead Measurement:
Place the probe in the middle of the eyebrows. Press the measurement button and slowly move the probe on the skin for about 3 seconds. Release the button when you hear the notification sound.
4. Measurement results will be displayed on the screen and recorded automatically.

Note:

When environmental temperature is low, forehead skin temperature exists rather big

difference, or probe fails to direct at eyebrows. This might cause the result that Ear measurement is lower than usual level.

To solve the problem, please press and hold the measurement button and move the probe between eyebrow (or temple). After having confirmed the probe scanning the eyebrow(or temple), release the button and finish the measurement.

MEMORY FUNCTION

Smart Thermometer can automatically store the latest 30 sets of measurement data of Ear and Forehead temperature.

Note: if you want to check recorded temperature data, please make sure the Smarte Thermometer is power off.

1. Press measurement button till the screen show "MEMORY" symbol before you check recorded temperature data.
2. Press measurement button continuously to view the records in turn.

POWER OFF

1. Press and hold the "power" button until the screen shutdown automatically.
2. This product will automatically shut down if there is no operation in one minute.

BLUETOOTH SETTING

1. Download Koogeek app: download link <http://app.koogeek.com/download>.
2. Install and log in: please register with an email account and log in.
3. Add a device: launch the app and enter the Home Page. Click "+" button to add the thermometer.
4. Click "install new device" and "Thermometer" button to connect.
5. Please follow the app instruction to finish remaining steps.

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4. Finalice la instalación de koogeek-T1.

Caution

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Manufactured in accordance with International Standards

Manufacturer

Refer to the instructions

Equipment of BF

Recycling and Processing

The european Union flag

Keep dry

The represent of European Union

EC-Representative :

Wellkang Ltd t/a Wellkang Tech Consulting Suite B, 29 Harley Street,
LONDON W1G 9QR, United Kingdom.

App Download Address

Contact us :

<http://www.koogeek.com/>

Email:support@koogeek.com